



PAL Handpiece | PAL-750

- 35% faster procedure¹
- 49% less surgeon fatigue¹
- Lasts 2X as long as PAL-650
- 2-year warranty
- 5 years of complimentary loaner equipment during repairs²
- Weight: 1.1 lb (0.5 kg)
- Stroke length: 3.2 mm
- Maximum speed: 4,000-5,000 cycles/min

1. Compared to traditional manual liposuction. Katz BE, Bruck MC, Coleman WP. The Benefits of Powered Liposuction Versus Traditional Liposuction: A Paired Comparison Analysis 2001. Dermatol Surg; 27:863-967.
2. For customers in the United States.

NEW! Redesigned PAL Cable

PAL Cable | 7020-PAL

- More lightweight and flexible³
- Thicker, more robust pins³
- Easier alignment and connection to the new PAL Console³

3. vs. 5006-PAL Cable



Aspiration Tubing | PAL-1200 Infiltration Tubing | PAL-INF-1600



Vacuum Tubing | ASP-TUBING




Biofilter | ASP-FLTR2



LipoFilter | ASP-CAN-2S

- 2,500 mL lipoaspirate volume
- Gravity-decanting fat processing system
- Sterile-packed, single-use

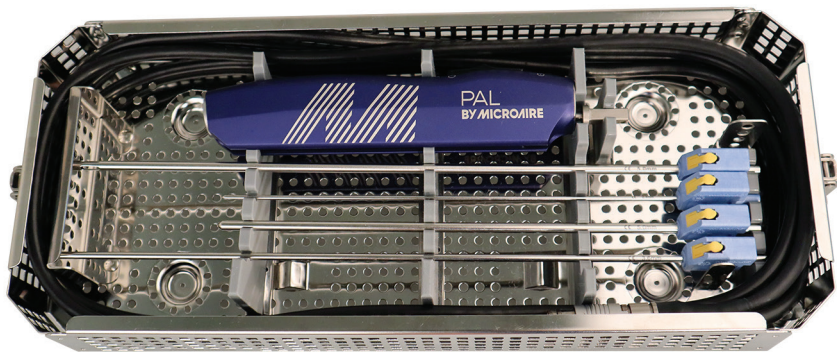

Sterilization Case | PAL-CASE

- Holds 1 PAL handpiece, 1 cable, 4 cannulas up to 30 cm in length
- Dimensions: 15.9" x 6.5" x 2.4"

To order call toll free **800.722.0822**



©2025 MicroAire Surgical Instruments, LLC
All rights reserved. LIT-PALSYNERGY Rev B

www.microaire.com

MicroAire Surgical Instruments
3590 Grand Forks Boulevard
Charlottesville, VA 22911 USA
800 722 0822 or +434 975 8000

PALSYNERGY®
More capabilities. More control.



MICROAIRE®

PAL SYNERGY®

More capabilities. More control.

PAL Synergy® gives you more control of your devices.

Run the infiltration and aspiration pumps using the multi-functional foot pedal or control them directly from your PAL Handpiece. Seamlessly transition between infiltration and aspiration without relying on additional hands.

PAL® Liposuction Canisters

PAL-SC-2000 | PAL-SC-3000

- Made with lipid-resistant materials
- 29 inHg (736 mmHg) vacuum pressure rating

Infiltration Pump

- Average max flow rate: 1,155 mL/min

Cart | CART

- Footprint: 28.2" x 27.2" (71.6 cm x 69.1 cm)
- Height (without poles): 39.2" (99.6 cm)
- Includes IV Pole, Canister Pole, Basket, 2 Canister Clamps and Rings

PAL® Console | 7020

Enables full system integration with connection points for:

- 2 PAL handpieces
- 2 multi-functional foot pedals
- Aspiration and infiltration pumps

Aspiration Pump

- Pump type: 1/3 horsepower, rocking piston
- Powerful vacuum pressure: >29 inHg (>736 mmHg)
- Quiet operation: 56.8 dB
- Customizable preset pressure button
- Internal fluid sensor and catch canister to prevent motor damage



System reverse view

- Color-coded connections for easy setup
- Single cord-to-wall plug powers the entire system⁴

4. When used with three power cord jumpers (ASP-CORD)

Multi-Functional Foot Pedal

- Select foot pedal mode and use the pedal to activate pumps
- Select PAL mode and use the pedal to select a pump for PAL handpiece activation

ORDER NOW!

SYNERGY-KIT

Includes:

- Infiltration Pump
- Aspiration Pump
- Multi-Functional Foot Pedal
- All required connection cables
- Infiltration Tubing
- 2 Biofilters

